



NEO BEAUTY BCN, S.L.U  
Oficina: AVDA. CARRILET, 243, PLANTA 1, PUERTA 3B,  
L'HOSPITALET DE LLOBREGAT (08907 BARCELONA)  
[www.neobeautybcn.com](http://www.neobeautybcn.com)  
[info@neobeautybcn.com](mailto:info@neobeautybcn.com)  
+34 629 523 166

## RADIO FREQUENCY MONOPOLAR



Art. IB-RF02

## NOTICE

**Attention Before operation, please read below information carefully.**

### Treatment attention:

1. Select proper treatment head according to specific body part.
2. The intensity should be adjusted from weak to strong, to avoid Sudden strong stimulation that makes the client get unnecessary fear
3. When make treatment, if the body feel pain or the skin become red and imprint, should stop the treatment at once
4. The treatment head should move slowly on the surface of skin, don't leave on any body part for too long time.
5. During treatment, it is forbidden to wear any metallic article.

### Taboo:

1. It is forbidden to apply treatment to pregnant women.
2. It is forbidden to apply treatment to epileptic and psychopath.
3. If people having fever, infection, skin disease(whelk,tetter,herpes) fresh scar wound, serious chronic disease(cancer,AIDS,hepatitis, blood disease etc.)or cruor disease, heart disease, it is for bidden to be applied the treatment.
4. If people having phlebitis, having medical operation in two months(vein cirrhosis),having varicosity disease or limb trouble, it is forbidden to be applied the treatment.
5. If people taking medicine such as anti inflamed medicine, anticoagulant, antibiotics and so on for a long time, and the skin is easy to be bruise after light impact, please inquiry the doctor before applying the treatment.


## FUNCTION AND PRINCIPLE

### ● Skin therapy

**Function:** rise and frap skin, whiten and tender skin, control dark Sore, reduce wrinkle, frap loose belly after giving birth.

**Principle:** Double pole radio frequency beauty function, it employ frequent swtich current. When the current pass through body tissue, it can make effect to the electron and generate the phenomena polarization electron movement and make electron impact each other, rotating and twisting, which can generate biology heat energy and increase the temperature deep inside of the skin. In this way to Promote collagen fibre to frap and reborn collagen albumen, therefore can remove wrinkle and frap and tender skin.

**Curing symptom:**

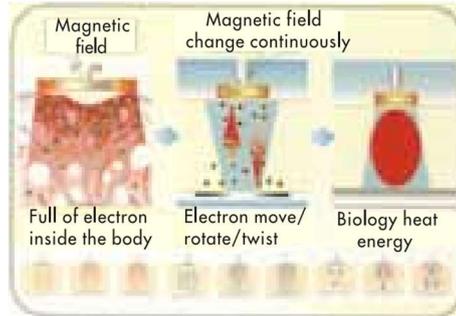

RF polarization working principle picture


1. Crow-feet
2. Blepharoptosis
3. Sagging of face and chin
4. Other sin ageing symptom

● **Fat therapy**

**Function:** Burning fat, slimming body, improve astriction, strengthen elasticity.

**Principle:** Double pole frequency resonate with body biology electricity, make fat cell break down, frap and release, meanwhile promote lymph system detoxin and accelerate blood circulation, improve metabolism.

**Curing symptom:**

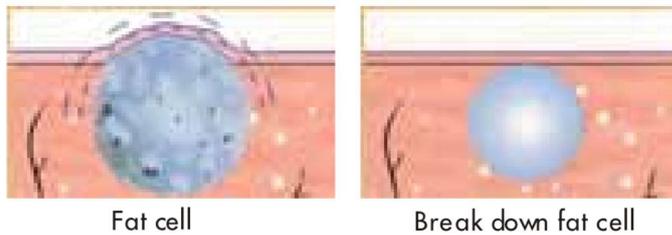


**TREATMENT EFFECT**




## Advantage:

|              |                       |
|--------------|-----------------------|
| No operation | No wound              |
| No intrude   | No need convalescence |
| No blood     | Short treatment time  |
| No pain      | Rapid effect          |
| No scar      | Long keep             |

## PANEL AND ACCESSORY

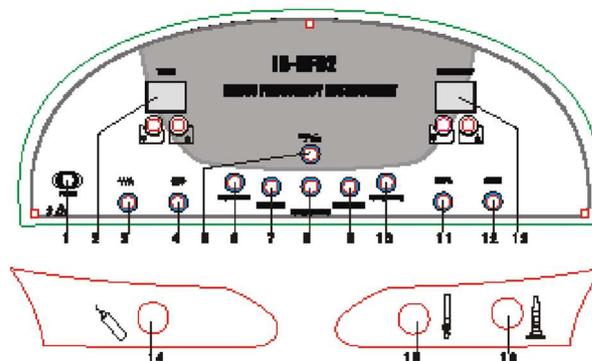

- |                           |                           |
|---------------------------|---------------------------|
| 1. Power                  | 9. Removes belly wrinkles |
| 2. Time display           | 10. Body slimming         |
| Time adjust               | 11. For face care         |
| 3. Continuous wave select | 12. For body care         |
| 4. Interval wave select   | 13. Intensity display     |
| 5. Start/pause            | Intensity adjust          |
| 6. Removes face wrinkles  | 14. Negative output port  |
| 7. Face lifting           | 15. Face care output port |
| 8. Detaches lifting       | 16. Body care output port |

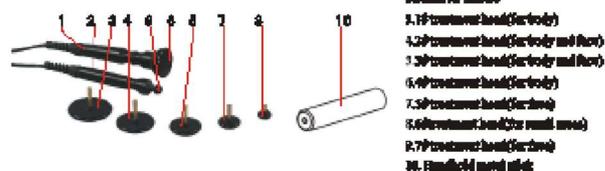

### Accessory

- 1. Big RF handle
- 2. Small RF handle
- 3. 18° treatment head (for body)
- 4. 30° treatment head (for body and face)
- 5. 30° treatment head (for body and face)
- 6. 40° treatment head (for body)
- 7. 50° treatment head (for face)
- 8. 60° treatment head (for small area)
- 9. 70° treatment head (for face)
- 10. Handhold metal stick

## OPERATION STEP AND ATTENTION

### Operation step:

- Clean skin.
- Connect the instrument to power supply.
- Connect the treatment heads to the instrument. Start the instrument.
- Client lies flat, and hold the metal stick.
- Spread beauty cream on the surface of treatment parts evenly.
- According to treatment part, choose proper treatment head, select the following output setting.

- Put the button FACE WRINKLE, time display 15 minutes, intensity display 11;
  - Put the button FACE LIFTING, time display 15 minutes, intensity display 9
  - Put the button BUTTOCKS LIFTING, time display 30 minutes, intensity display 8;
  - Put the button BELLY WRINKLE, time display 45 minutes, intensity display 10
  - Put the button BODY SLIMMING, time display 45 minutes, intensity display 8.
7. Put the button START/PAUSE, the instrument start working, time decrease gradually, and The blue sall point on the screen begin to flash.
  8. After choose the output mode, the working time and intensity can also be adjusted according to specific demand.
  9. Time goes to zero, the operation is end.
  10. Clean and sterilize the accessories.

### **Treatment attention:**

1. During operation, because of individual difference, so need to communicate with client and know their feeling ,should adjust the intensity according to specific demand.
2. For the treatment of sensitive part, can adjust the intensity to 3 to 5,which is subject to client's specific condition.
3. For the small treatment head ,the current is concentrated, the intensity should be adjusted to smaller; For the big treatment head, the current is dispersed, the intensity should be adjusted to bigger
4. Treatment head should conitnuous move upon skin, should not stay In a part for more than 3 seconds.
5. If continuous using for half an hour, should stop the operation and rest it for at least 10 minutes for the next operation, in this way can extend the life of Instrument.

## OPERATION GUIDE

**Attention: During treatment, there must be no direct body touch between the client and any other people, the beautician should wear insulating glove, in this way to avoid unnecessary injury.**

1. Picture guide for face

A. Lift chin, improve mouth's feet  
(Use 2# 3# 4# 5# treatment head)

B. Lift cheek, improve face loose  
(Use 2# 3# 4# 5# treatment head)

C. Improve crow's feet  
(Use 3# 4# 5# treatment head)

D. Improve forehead wrinkle  
(Use 2# 3# treatment head)

E. Get rid of eye pouch and dark eye socket.  
(Use 6# 7# treatment head)

2. Picture guide for belly

(Use 1# treatment head)

F. Clockwise movement

G. Anticlockwise movement

H. Movement from down to up

3. Picture guide for thigh

(Use 1# treatment head)

I. Movement from inside to outside, from down to up

4. Picture guide for breast lifting


(Use 1# treatment head)

K. Lift from down to up

5. Picture guide for hips lifting

(Use 1# treatment head)

J. Lift from outside to inside, from down to up

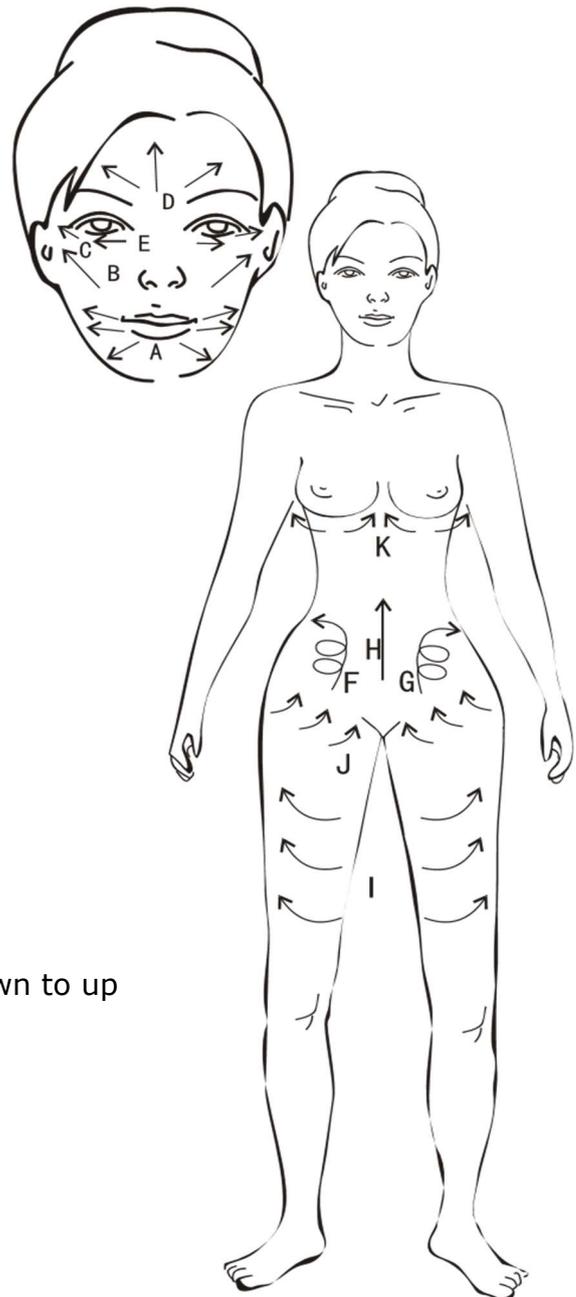

## TECHNICAL PARAMETER AND SOLUTIONS FOR SIMPLE FAILURE

| FAILURE   | SOLUTION   |
|---|--|
| <b>POWER ON, THE MACHINE IS NO RESPONSE</b>         | check whether the power supply is in good connection with main machine                           |
|   | check whether the power switch is on   |
| <b>NO OUTPUT</b>                                    | check whether the accessory is in good and right connection with the output port of main machine |
| <b>DURING WORKING, THE MAIN MACHINE IS OVERHEAT</b> | using For too long time, Power off and naturally cool Down                                       |
|   | check whether electricity voltage is correct (need Professional People to check)                 |
| <b>NO POWER</b>                                     | check the fuse, if broken need to replace it.  |

| VOLTAGE   | POWER | HZ.      | SIZE            |
|-----------|-------|----------|-----------------|
| 240/110 V | 100 W | 50/60 Hz | 54 x 40 x 23 cm |

## PACKING LIST

| NAME            | UNIT  | QTY | REMARK            |
|-----------------|-------|-----|-------------------|
| Main machine    | set   | 1   |                   |
| Big RF handle   | Pc    | 1   |                   |
| Small RF handle | Pc    | 1   |                   |
| Metal stick     | Pc    | 1   |                   |
| Treatment head  | set   | 1   | 7 treatment heads |
| Power wire      | Strip | 1   |                   |
| Fuse            | Pc    | 2   |                   |
| Instruction     | Pc    | 1   |                   |

## NOTICE

1. The instrument has been delivered with full items, when received the instrument, please check all the accessories whether something is missing, and refer to our packing list.
2. If something is broken during transportation, please report to the dealer and send us a copy information.
3. Please reserve the original packing items, when need maintain, it is necessary to send us the instrument with original packing.

## SAFETY NOTICE

1. The machine is professional instrument for beauty salon, in order to avoid any unnecessary hurt, need professional beautician to make operation.
2. Do not disassemble the machine or operate the machine that not explained in the instruction book. All the maintenance services shall be carried out by the professional maintenance technician approved by our company.
3. Do not install or operate the machine with wet hands or near the water sources so as to avoid the liquid spraying/swashing over the machine.
4. Pull out the plug and connect with the maintenance service center in the following cases:
  - A. liquid spray/swash over the machine or the machine get caught in rain or water.
  - B. the machine has peculiar smell, fume/ smog or abnormal noise.
  - C. power wires damaged or scratched.
  - D. instrument fell down or the shell damaged.
5. Do not put any item on the power wires, or install the instrument in the place where its power wires will be trampled.
6. Pull out the power wires at the end of use for the sake of safety.
7. Do not insert any type of item into the instrument opening; otherwise there will be a fire or electric shock. In case of something fell into the opening but you could not take it out, please pull out the power wires and connect with the dealer or our company.
8. Do not put the instrument on the unstable trolley, shelf or vehicle, if fall down, the instrument may be damaged.
9. To guarantee the life circle of this machine, the interval of the continuous switch on and off should be over one minute.
10. The power shall comply with the specification specified in the label, otherwise the machine will not run or run abnormally, even be burnt out.
11. The company has right to revise or amend this manual.

The company reserves all the rights of interpreting the information mentioned above.

**Attention:**

If the supply cord is damaged, it must be replaced by manufacturer, its service agent or similarly qualified persons in order to avoid a hazard.

This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.

Children should be supervised to ensure that they do not play with the appliance.